

CleanCut® & PerfectCut®

AORTIC PUNCHES

For many surgeons, the aortic punch is the perfect anastomotic device. We offer surgeons a wide variety of aortic punches, including a unique aortotomy system that is ideal for creation of uniform openings at the aortotomy site for vein grafts, it's called the PerfectCut®. The PerfectCut® system is designed to create a precise, consistent, reliable, repeatable incision and circular hole that eliminates lateral side notches from the aortotomy. The PerfectCut® and CleanCut® punches create consistent openings at the anastomotic site through a rotating, self-aligning cutting mechanism that produces uniform openings through precise cutting actions with reduced force.

FEATURES & BENEFITS

- Produces uniform openings at the anastomotic site
- Rotating cutting action requires less force
- Self-aligning cutter mechanism
- Six sizes, two lengths, and three anvil shape options: bullet-nose, standard, and flat

SPECIFICATIONS

- Materials:
 - Medical-grade stainless steel
 - Thermoplastic construction
- Available sterile and non-sterile directly from our manufacturing facility and from all procedure kit suppliers

REF	Size	
RCL28, RCC28	2.8 mm	●
RCL35, RCC35	3.5 mm	●
RCL40, RCC40, RCB40	4.0 mm	●
RCL45, RCC45, RCB45	4.5 mm	●
RCL48, RCC48, RCB48	4.8 mm	●
RCL50, RCC50, RCB50	5.0 mm	●
ALX33 (Lancet)	3.3 mm	+
<i>Punches 6 units/case Lancet 25 units/case</i>		

Manufactured by Quest Medical, Inc.

One Allentown Parkway | Allen, Texas 75002-4211 USA 
 custserv@nordsonmedical.com | +1.800.627.0226 | +1.972.390.9800

© 2025 Nordson MEDICAL. All rights reserved. PerfectCut® is a registered trademark of Nordson MEDICAL.

For complete indications, contraindications, warnings, precautions, and potential adverse events, refer to the product's Instructions for Use (IFU).

ADVANCED TOOL SET FOR PERFECT AORTIC INCISIONS

PERFECTCUT® AORTIC LANCET

- Consistent, reliable, repeatable, crisscross incisions
- For aortotomies such as at the coronary artery proximal anastomotic and arterial cannulation sites
- Lancet blade features retractable safety shield
- 6.5 inches (16.51 cm) with 10 x 3.3 mm blade


PERFECTCUT® ROTATING BULLET-TIP PUNCH

- Punch utilizes the lancet's uniform incision to reduce the likelihood of lateral side notches caused by inconsistent scalpel blade incisions
- Advanced "Bullet-Tip" designed for straight insertions— no "button-hole" technique necessary
- Rotating cutting action requires less force and eliminates shearing of delicate tissue
- Self-aligning cutter mechanism minimizes mechanical jamming
- 8 inches (20 cm)

REF	Size
RCK40	4.0 mm + ●
RCK45	4.5 mm + ●
RCK48	4.8 mm + ●
RCK50	5.0 mm + ●
Kits sold 6 units/case	

CREATING THE PERFECT AORTOTOMY WITH THE PERFECTCUT® SYSTEM

- Step One:**
Aortic lancet creates precise crisscross incisions
 - Step Two:**
Bullet-Tip Punch inserts easily into incisions
- ➔ **Result:**
The solid, durable construction results in reliable performance and a perfect aortotomy

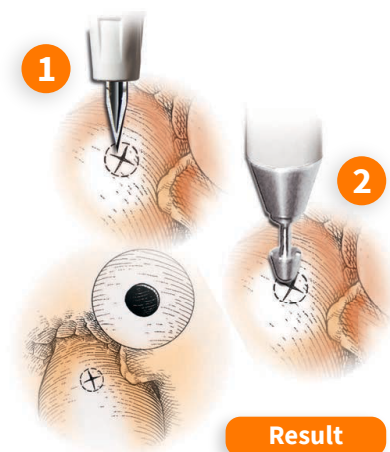

The CleanCut® and PerfectCut® Aortic Punches are intended for use in cardiac surgery to create aortotomies for saphenous vein graft anastomosis during coronary artery bypass graft (CABG) procedures.