

Retract-O-Tape® & SaddleLoop®

VASCULAR LOOPS

Retract-O-Tape® and SaddleLoop® are sealed silicone tubes that provide an air cushion for occluding or retracting vessels during surgery, with minimum trauma. All varieties maintain the round profile while providing necessary flexibility and are latex-free.

Our vascular loops are available both sterile and non-sterile directly from our manufacturing facility and all procedure kit suppliers.

RETRACT-O-TAPE®

- Occludes or retracts a vessel without risk of tissue damage
- Maintains round profile while providing necessary flexibility
- Sealed silicone tubing provides an air-cushion
- Needle version offers superior control for safe, quick positioning

SPECIFICATIONS

- Materials:
 - Needle: surgical stainless steel
 - Tubing: radiopaque silicone
 - Latex free
- Available sterile and non-sterile from Nordson MEDICAL and from all procedure kit suppliers

Manufactured by Quest Medical, Inc.
 One Allentown Parkway | Allen, Texas 75002-4211 USA
 custserv@nordsonmedical.com | +1.800.627.0226 | +1.972.390.9800
 © 2025 Nordson MEDICAL. All rights reserved. Retract-O-Tape® and SaddleLoop® are registered trademarks of Nordson MEDICAL.

REF

1041	Blunt needle vascular loop, specially designed, pre-attached, curved blunt tip 18G needle; thin wall Length: 12 in. (30.48 cm) Diameter: 0.049 in. (1.24 mm)
1015	Retract-O-Tape® vascular loop; thin wall Length: 18 in. (45.72 cm) Diameter: 0.131 in. (3.32 mm)
1012	Retract-O-Tape® vascular loop; thin wall Length: 18 in. (45.72 cm) Diameter: 0.066 in. (1.67 mm)
1032	Retract-O-Tape® vascular loop; thick wall Length: 18 in. (45.72 cm) Diameter: 0.066 in. (1.67 mm)
1001	Retract-O-Tape® vascular loop; thin wall Length: 12 in. (30.48 cm) Diameter: 0.049 in. (1.24 mm)
1002	Retract-O-Tape® vascular loop; thin wall Length: 12 in. (30.48 cm) Diameter: 0.066 in. (1.67 mm)
1071	Semi-pointed needle vascular loop, specially designed, pre-attached, curved blunt tip 18G needle; thin wall Length: 12 in. (30.48 cm) Diameter: 0.049 in. (1.24 mm)

*Sold in 10 packets of 2 units (total: 20 units/case) except
 REF 1015, 10 packet of 1 unit (total: 10 units/case)

Retract-O-Tape® & SaddleLoop®

VASCULAR LOOPS



BETTER VISIBILITY WITH MINIMAL INTRUSION

SADDLELOOP®

- Radiopaque silicone tubing with sealed ends provide an air-cushioned surgical tool for manipulating anatomical organs and structures or ligating vessels with minimal trauma
- Materials:
 - Needle: surgical stainless steel
 - Tubing: radiopaque silicone
 - Locking mechanism: ABS with barium sulfate
 - Latex free
- Available sterile and non-sterile directly from our manufacturing facility and from all procedure kit suppliers

REF


1541 Blunt, 18G needle tip | Length: 4 in. (10.16 cm)

1571 Semi-pointed, 18G needle tip | Length: 4 in. (10.16 cm)

Sold in 8 packets of 2 units (total: 16 units/case)

SADDLELOOP® TECHNIQUE

The SaddleLoop® device is perfect for minimally invasive procedures where manipulation is restricted and can be used within the confines of a heart stabilizer.

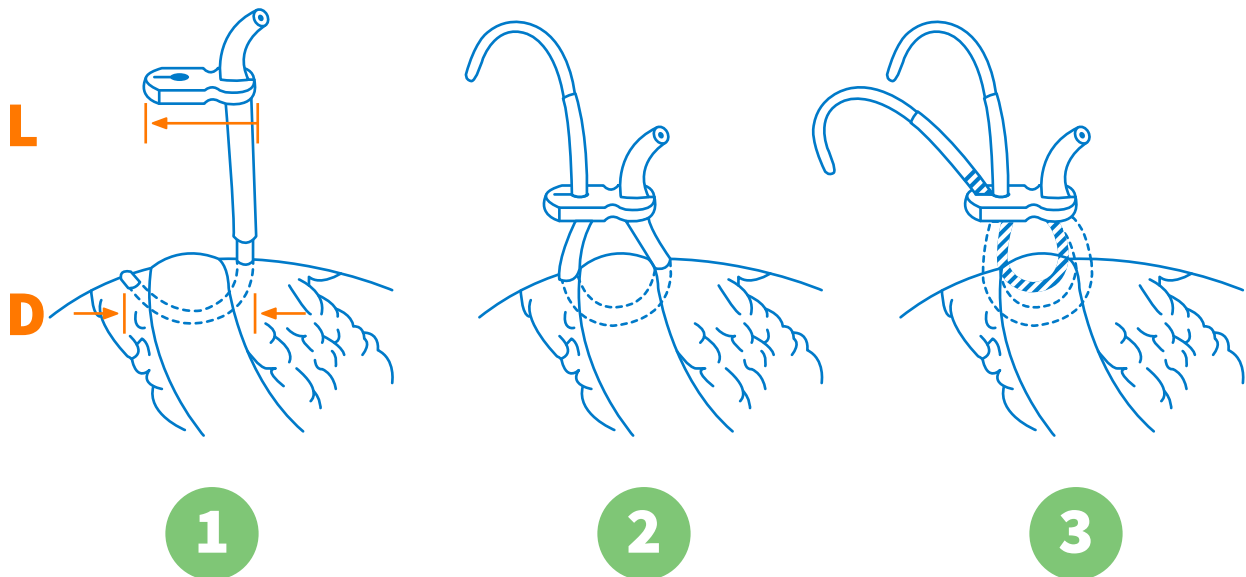

Ensure distance "D" is greater than length "L"

Push needle through locking mechanism

Slightly stretch tubing to tighten before sliding into slot of keyhole locking mechanism